


## 12 Fr x 45 cm

NOW AVAILABLE for use with the GORE® EXCLUDER® Iliac Branch Endoprosthesis






Hydrophilic coating and enhanced flexibility provide exceptional access to challenging anatomies and branch vessels.

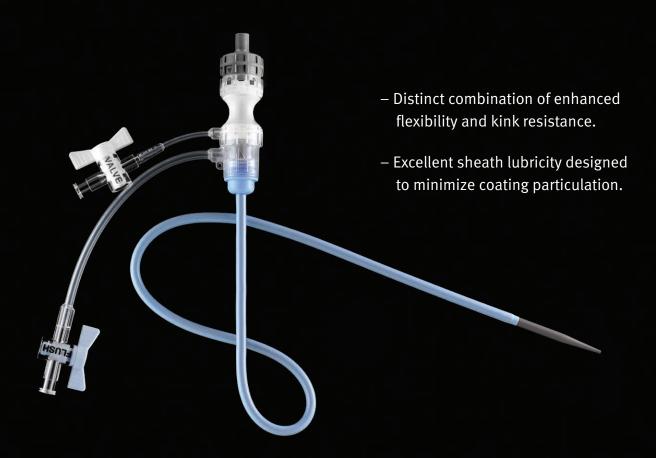



## Exclusive DrySeal valve enables introduction of multiple devices with proven hemostasis control.

- Accommodate multiple wires and catheters through the valve simultaneously with minimal blood loss.
- Facilitate hemostasis without intraprocedural manipulation of the valve.

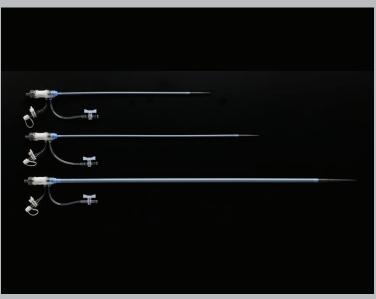


## Optimized profile and new configurations provide tailored delivery options for a broader range of patient anatomy.

- Decreased sheath wall thickness for larger inner diameter and decreased outer diameter.
- Additional working lengths of 33, 45,
  and 65 cm to care for more patients.



## Engineered for use with our endovascular portfolio.

 12 Fr x 45 cm configuration designed for and tested specifically with the GORE® EXCLUDER® Iliac Branch Endoprosthesis. The GORE® DrySeal Flex Introducer Sheath with hydrophilic coating consists of an introducer sheath with DrySeal valve attached, a dilator, and a 2.5 ml valve inflation syringe.

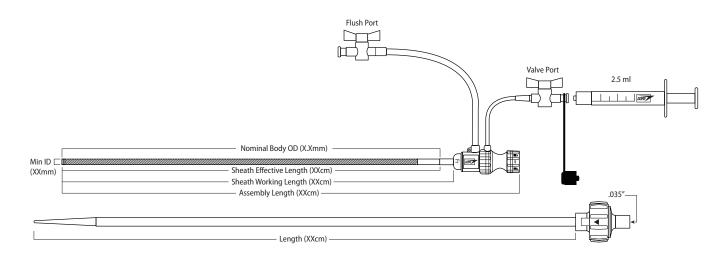

	Configu	IRATION	Мінімим	Nominal	SHEATH Effective	Sheath Working	Assembly	DILATOR
Catalogue Number	SHEATH SIZE (Fr)	LENGTH (cm)	SHEATH ID (mm)	Body OD (mm)	LENGTH (cm)	LENGTH (cm)	LENGTH (cm)	LENGTH (cm)
DSF1233	12	33	4.0	4.7	32	33	39	42.9
DSF1245	12	45	4.0	4.7	44	45	51	55.1
DSF1433	14	33	4.7	5.3	32	33	39	43.7
DSF1533	15	33	5.0	5.6	32	33	39	43.9
DSF1633	16	33	5.3	6.1	32	33	39	44.2
DSF1833	18	33	6.0	6.7	32	33	39	45.0
DSF2033	20	33	6.7	7.5	32	33	39	45.5
DSF2065	20	65	6.7	7.5	64	65	71	77.5
DSF2233	22	33	7.3	8.2	32	33	39	46.2
DSF2265	22	65	7.3	8.2	64	65	71	78.0
DSF2433	24	33	8.0	8.8	32	33	39	46.7
DSF2465	24	65	8.0	8.8	64	65	71	78.7
DSF2633	26	33	8.7	9.5	32	33	39	47.5
DSF2665	26	65	8.7	9.5	64	65	71	79.5



W. L. GORE & ASSOCIATES, INC. Flagstaff, AZ 86004

+65.67332882 (Asia Pacific) 00800.6334.4673 (Europe) 800.437.8181 (United States) 928.779.2771 (United States)

goremedical.com

Consult Instructions for Use

Products listed may not be available in all markets.